

Donor Doc *Phlebotomy*TM

Electronic Blood Draw Data



Improve Compliance.
Enhance Efficiencies.
Help Reduce Errors.

New Features and Benefits:

Safety


- Presents the donor's CMV status on the Donor Verify and Donation Summary pages
- Captures and records a donor's reaction to the phlebotomy process, selects treatment, displays physical exam values and records comments

Flexibility

- Reconnects to the main server and uploads information collected from an offline configuration used at mobile drives
- Offers optional pages to capture saline and anticoagulant used during the draw for procedures that require use of those supplies


Efficiency

- Records red cell and plasma loss values in milliliters within a new page of the workflow
- Allows the user to define collection volume per gender, weight and collection type by searching the target volume, and calculating the ideal draw for the donor

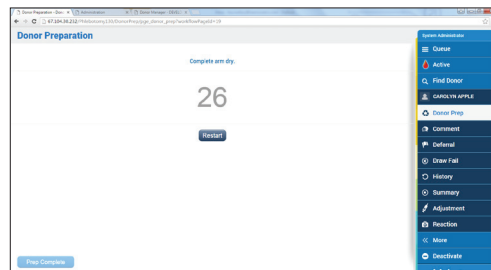

Configures two timers while preparing the arm for phlebotomy, including scrub time and dry time within the Time Management Settings Table.


Supports Dynamic Wild Card, allowing the user to customize the label and capture supplies and equipment used in the phlebotomy process.

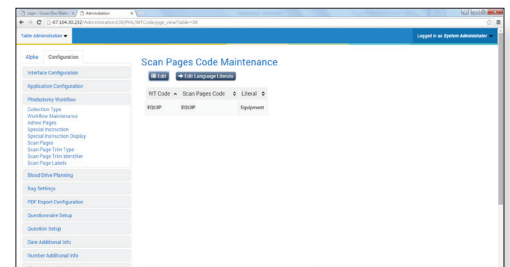

Donor Doc Phlebotomy™ software enables greater comfort, accuracy, flexibility and reliability to the phlebotomist through a paperless, compliant donation process. The flexible user-defined workflow streamlines data, and provides real-time active donor status. By eliminating manual entry, maintenance and paper records, the software helps to improve the donation experience for both the donor and the user.


Phlebotomy Queue

- Displays a list of all pre-qualified donors ready to be drawn
- Navigates to the appropriate forms and verification page for the next selected donor
- Updates the donor collection type automatically if modified prior to the draw

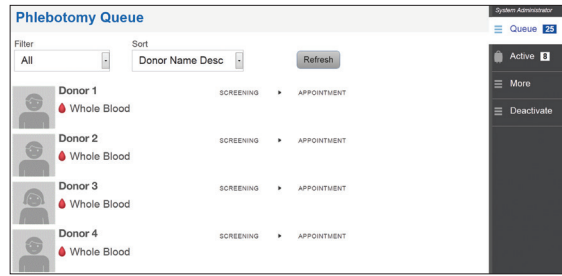


Active Queue

- Provides a list of donors who have already been selected from the phlebotomy queue
- Allows user to switch back-and-forth between donors they are actively working with
- Displays the donor's current phlebotomy step
- Utilizes a donor graphic so the user can select the arm(s) used for phlebotomy within the Arm Selection page

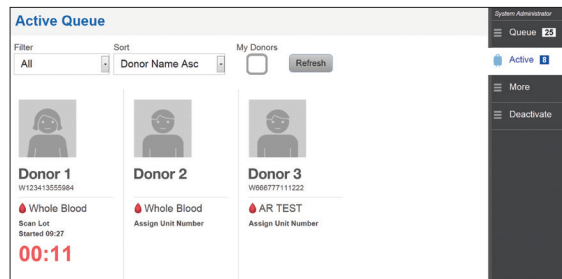


End Donation

- Allows user to determine if the draw was successful or indicate failure code. Deferral codes can also be defined for any failed code
- Captures issues related to blood draws
- Records actual volume collected for each main bag during a draw, with the ability to compare actual to intended volume

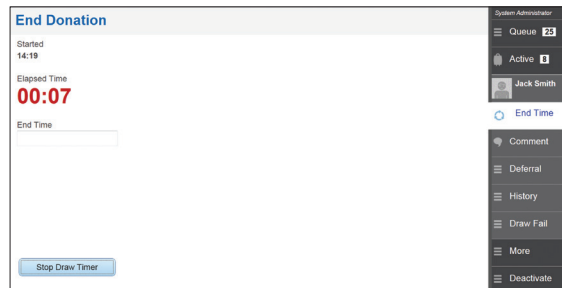

Donation Summary

- Displays donor's phlebotomy summary and allows users to complete phlebotomy and submit data to host system
- Allows User, at a glance, to edit data if necessary before sending to host
- Prints manufacturing label with user configurable information following the release of the record

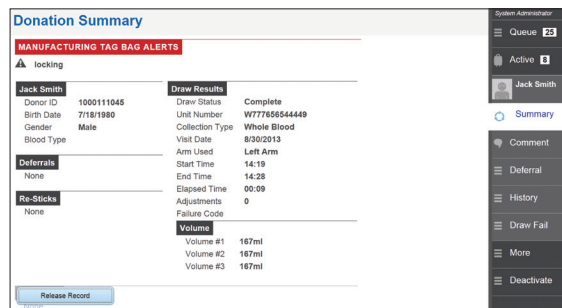

To learn more, visit our website at www.haemonetics.com or contact your local Haemonetics Account Manager.

For a list of worldwide office locations and contact information, visit www.haemonetics.com/officelocations

© 2016 Haemonetics Corporation. Haemonetics and Donor Doc Phlebotomy are trademarks or registered trademarks of Haemonetics Corporation in the USA, other countries, or both. All rights reserved. 06.2016 USA. COL-COPY-000743(AC)